


Light therapy Lamp

User Manual

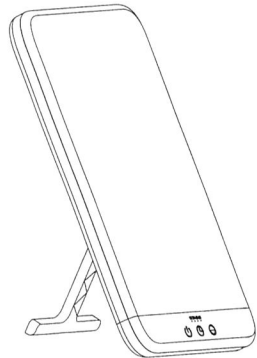


Introduction

Why use a light therapy lamp?

Natural spectrum light therapy lamp may help to reduce the effects associated with winter blue, jet lag, shift work and season time changes and boost your sense of well-being.

Product diagram

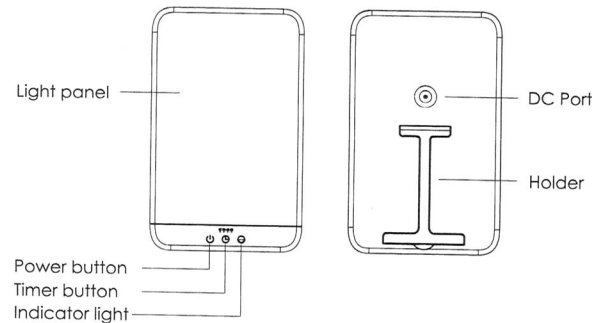

Specifications

Power	12W
Working voltage	DC12V, 1.2A
Lighting source	LED SMD2835 0.5W 96Pcs
CCT	6000K/2300K
CRI	Ra>92
Light intensity	10000lux@10cm
Materials	ABS+PC (Fireproof)
Dimensions	L240*W150mm*H16.5mm

How to use

Use the power adapter to connect the lamp with wall outlet.

1. Power ON/OFF & Brightness control

Tap the power button to turn on the lamp at 100% brightness level. Tap the power button in turn to adjust another three brightness levels (75%, 50%, 25% brightness level). Tap again to turn off the lamp. When the lamp is ON at any brightness level, long press the power button to turn it off. Only in this case the lamp will remember the last setting before turned off.

2. 15/30/45/60min Auto-timer

Tap the timer button in turn to activate the 15/30/45/60minutes auto-timer, the white indicator light will be constantly on. Tap it again to cancel the timer.

Package Contents

- 1x Light therapy lamp
- 1x Power adapter
- 1x User manual

Warning

1. The light therapy lamp is only intended for use on the human body.
2. Connect the device only to the mains voltage listed on the type plate.
3. Do not dip the device into water and do not use it in wet rooms.
4. Do not use it in the vicinity of people who are helpless, children under age 3 or people who are insensitive to heat (people with skin lesions due to illness); do not use after taking painkillers or drugs or drinking alcohol.
5. Do not touch the device with wet hands while it is plugged in; do not allow any water to be sprayed onto the device. The device must be operated only when it is completely dry.

6. Do not use the device if it shows signs of damage or dose not function properly. In these cases, contact our customer services.
7. You can not repair the device. It contains no parts that you can repair.
8. Do not make any changes to the device without the manufacturer's permission.

Warranty & Customer Support

This lamp is covered with 12 months warranty from the date of its original purchase from authorized sellers. For questions, support, or warranty claims, contact our customer services for help.

Statement

We can only provide after-services for the products that are sold by authorized sellers. if you have purchased from an unauthorized seller, please contact them directly for service or warranty issues.



The touch button area will appear heat phenomenon when the lamp was used for a long time. This is a normal phenomenon, it will not affect the normal function.